

CHEMGLASS

Air-Operated Drain Valves for 10 through 50L Chemglass Reactors

- Fits existing Chemglass reactors with 2" beaded pipe drains
- Ideal for solvent extractions & partitions after a reaction
- Great for low temp reactions where standard valves need constant tightening to prevent leakage
- Held closed by spring force & opened by air pressure
- Chemglass Part Nos. **CG-1968-V-series**



Control Panel

Included are two, factory-preset regulators. The single control knob quickly opens and closes the valve as needed. Set-up is easy and flexible as the panel may be mounted on the stand or outside a fume hood.

Please Note: Due to the automatic nature of this valve, we do not recommend it for slurries or viscous materials.



Bottom Valve

Features a 3/4" compression fitting on the drain outlet. Nose and housing are PTFE and the shaft is available in 316L stainless or Hastelloy®. The valve mounts easily on existing 2" beaded pipe drain.



Optional Secondary Valve Kit

For emergency shut-off of the reactor contents in the unlikely event of valve failure.



See back page
for ordering
information. ►

Air-Operated Drain Valves for 10 through 50L Chemglass Reactors




Drain Valves with Stainless Steel or Hasteloy® Shafts

| Part Number | Description | Price |
|--------------|--|--------------------|
| CG-1968-V-01 | Air-Operated Lower Drain Valve, Stainless Steel Shaft , Complete with Air Control Panel, Perfluoro O-rings, Tubing and Fittings | *\$3,085.45 |
| CG-1968-V-02 | Air-Operated Lower Drain Valve, Hasteloy® Shaft , Complete with Air Control Panel, Perfluoro O-rings, Tubing and Fittings | *\$4,131.40 |

Mounting clamp not included. See **CG-9253** Universal Kwik Clamp.

Please Note: Due to the automatic nature of this valve, we do not recommend it for slurries or viscous materials.

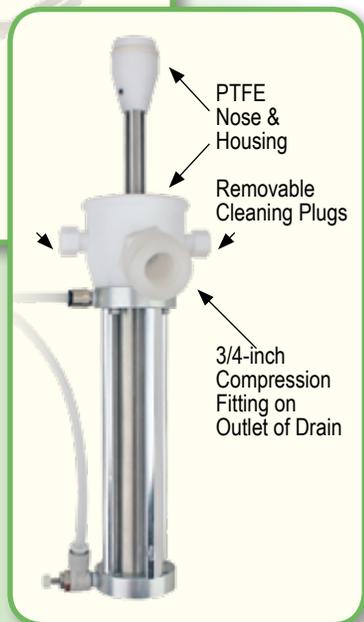

See it in action on VIMEO.com

Search for: Chemglass

It's also on the product page of our website.

Search for: CG-1968-V



Optional Secondary Valve Kit

| Part Number | Description | Price |
|--------------|--|--------------------|
| CG-1968-V-10 | Secondary Valve Kit, All-PTFE Ball Valve with 3/4" Compression Fittings for Use with Air-Operated Drain Valves | †\$1,474.75 |



Hasteloy is a reg. trademark of Haynes International, Inc.