

HandyStep[®] touch

The next generation repetitive pipette

Simple

Touchscreen operation

Easy

Automatic tip ejection

Efficient

Compatible with various third party tips



BRAND

BRAND. For lab. For life.



BRANDTECH[®]
SCIENTIFIC, INC.

To Place Orders: Call -800-843-1794 • Email: cs@cglifesciences.com • Visit: www.chemglass.com

HandyStep® touch · HandyStep® touch S

Overview

- Touchscreen operation with intuitive menu navigation and integrated help function
- Two models available: HandyStep® touch or HandyStep® touch S with expanded functionality
- Menu languages: English, German, Spanish, French and Chinese
- Variable volume adjustment from 1.0 μ L to 50 mL
- Inductive charging with optional charging stand

A closer look

Versatile

HandyStep® touch with the standard modes Multi-Dispensing, Auto-Dispensing, Pipetting and Favorites;

HandyStep® touch S with the additional functions Sequential Dispensing, Multi-Aspiration and Titration

Long-lasting

Many hours of continuous operation with a single battery charge

Robust

High-quality materials provide long service life, and smooth surfaces make cleaning/ disinfection easy

USB-C connection

For quick charging by cable

STEP button

Ergonomically arranged STEP button for effortless volume filling and delivery

Easy

Automatic size detection of BRAND® PD-Tip™ II

Efficient

In addition to BRAND® PD-Tip™ II, compatible third party tips can also be used.

Simple

Automatic tip ejection by the press of a button prevents contamination

Convenient

Inductive charging with the optional charging stand



Intuitive operation

Be flexible and have the right solution at hand for every application. Even complex procedures are easy to carry out – no problem with the expanded functionality of the HandyStep® touch *S*.

Mode	Function	HandyStep® touch	HandyStep® touch <i>S</i>
Multi Dispensing (MULTI-DISP)	Successive dispensing of a set volume	✓	✓
Auto Dispensing (AUTO-DISP)	Automatic dispensing by the instrument	✓	✓
Pipetting (PIP)	Pipetting a defined volume	✓	✓
Sequential Dispensing (SEQ-DISP)	Successive dispensing of different volumes	–	✓
Multi Aspirating (MULTI-ASP)	Aspirating different volumes	–	✓
Titrating	Manual or successive delivery of liquids for titration	–	✓
Favorites	Saving favorites for quick access	✓	✓
Settings	Individual customization of the device (e.g. brightness or device name)	✓	✓




The clearly structured menus make fast work with intuitive operation possible, while the most important functions and parameters are always in view.




Easy to understand
Main menu with all modes for quick selection



Clearly structured
Working mode with all important information

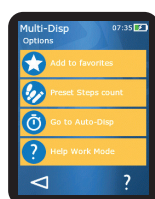

Organized
All options visible at a glance

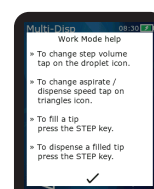

Helpful
Internal help right where it is needed

The perfect team

Automatic size detection with PD-Tip™ // from BRAND®

The HandyStep® touch saves time and prevents errors through automatic tip size recognition of the PD-Tip™ // from BRAND®, with patented coding on the pistons. After inserting the tip, the size is automatically recognized and displayed, making it easy to select the volume to be dispensed. When a new PD-Tip™ // of the same size is inserted, all instrument settings are maintained.



Multi dispensing with the HandyStep® touch and PD-Tip™ //



Sequential dispensing with the HandyStep® touch *S* and PD-Tip™ //

Technical data

HandyStep® touch HandyStep® touch S

Repetitive pipette, includes universal AC adapter (for USA, Europe, UK, Japan, Australia and New Zealand), shelf/rack mount, sample PD-Tip™ II, operating manual and performance certificate.

Description	BrandTech® Cat. No.	Chemglass Cat. No.	2021 List Price
HandyStep® touch	705200	CLS-3468-200	\$1,055.00
HandyStep® touch S	705210	CLS-3468-210	1,299.00

Inductive charging stand bundles

HandyStep® touch & Charging stand	705201	CLS-3468-200P	\$1,148.00
HandyStep® touch S & Charging stand	705211	CLS-3468-210P	1,525.00

Accessories

Description	BrandTech® Cat. No.	Chemglass Cat. No.	2021 List Price
Charging stand	705220	CLS-3468-300	\$351.80
Support stand	705230	CLS-3468-302	146.80
Shelf/rack mount	705235	CLS-3468-306	35.60
Universal AC adapter	705250	CLS-3468-304	86.60
Spare battery	705225	CLS-3468-308	81.60



PD-Tip™ II

Precision dispenser tips coded for automatic size recognition.

Capacity mL	Pack of	BrandTech® non-sterile Cat. No.	Chemglass non-sterile Cat. No.	2021 List Price	BrandTech® BIO-CERT® Cat. No.	Chemglass BIO-CERT® Cat. No.	2021 List Price
0.1	100	705700	CLS-3468-N-000	\$125.20	705730	CLS-3468-A-030	\$171.60
0.5	100	705702	CLS-3468-N-002	125.20	705732	CLS-3468-A-032	171.60
1.0	100	705704	CLS-3468-N-004	125.20	705734	CLS-3468-A-034	171.60
1.25	100	705706	CLS-3468-N-006	125.20	705736	CLS-3468-A-036	171.60
2.5	100	705708	CLS-3468-N-008	125.20	705738	CLS-3468-A-038	171.60
5	100	705710	CLS-3468-N-010	125.20	705740	CLS-3468-A-040	171.60
10	100	705712	CLS-3468-N-012	125.20	705742	CLS-3468-A-042	171.60
12.5	100	705714	CLS-3468-N-014	125.20	705744	CLS-3468-A-044	171.60
25	50 + 1 adapter	705716	CLS-3468-N-016	84.80	705746*	CLS-3468-A-046	52.60
50	25 + 1 adapter	705718	CLS-3468-N-018	39.80	705748	CLS-3468-A-048	52.60

Set of PD-Tip™ II, non-sterile, 20 pieces per size:
0.5/1/1.25/2.5/5/10/12.5 mL

705720 CLS-3468-N-050 169.60

Adapter for 25 and 50 mL
PD-Tip™ & PD-Tip™ II, PP,
autoclavable

Pack of 10

702398 CLS-3468-K-020 99.60

Pack of 5

702399 CLS-3468-K-010 78.40

*BIO-CERT® : Pack of 25 + 1 adapter



For more information visit www.brandtech.com or scan the QR Code.

HandyStep® touch, HandyStep® touch S and PD-Tips™ II have patents & patent applied for features.

California Residents: For more information concerning California Proposition 65, please refer to www.brand.de/calprop65.

Subject to technical modification without notice.



(800)-843-1794 | chemglass.com

888-522-2726 | www.brandtech.com

To Place Orders: Call -800-843-1794 • Email: cs@cglifesciences.com • Visit: www.chemglass.com